

Designed to protect



Product Features

- 1 LINE and LABEL**
Aligning the fill line on the lid with the label on the container ensures the internal safety flap inside the lid operates efficiently.
- 2 Simple FOUR CLICK assembly**
An audible click confirms the lid and base are assembled easily, correctly and securely.
- 3 Translucent lids**
A unique safety feature which indicates when the contents reach the fill line, without having to open or look into the aperture.
- 4 Clear segregation SIGNPOSTING**
✓ Waste type printed on the lid to assist with easier and safer segregation of sharps waste.
- 5 Large aperture with internal safety flap**
The large aperture with keyhole features ensures safe and easy disposal of sharps.
- 6 Temporary closure**
For increased user and patient safety when the container is not in use.
- 7 Fill line**
The inkjet printed fill line is clearly visible externally and in combination with the translucent lid and large 'headroom' ensures Sharpsafe® is the safest sharps disposal system available.
- 8 In-mould label**
This unique label is moulded directly into the container providing a permanent audit trail which cannot be removed.



New to the Sharpsafe® Range

The 25 L Sharpsafe®XL has been designed to contain laparoscopic instruments in the operating theatre. The lid has a fluid tight adhesive seal to prevent any leakage after final closure.

Features

- Fluid leak resistant seal with temporary closure feature
- Wall or trolley mountable for extra stability
- Available with a protected access insert
- Carry handle
- Interstackable for efficient space utilisation



Introducing Safe-Point®

The dental and podiatry needle remover.

Safe-Point® automatically unscrews and disposes of a used needle in seconds providing a responsible and effective solution to the safe disposal of used needles.

Safe-Point®		
SAFE-POINT UNIT For OCTAGONAL Needles	Cat No.	34300001
SAFE-POINT UNIT For ROUND Needles	Cat No.	34301001
SAFE-POINT UNIT For MORITA Needles	Cat No.	34302001
Each Safe-Point unit is supplied with 3 sharps containers.		
SAFE-POINT CONTAINERS 6 Replacement Sharps Containers	Cat No.	40002430



Sharpsafe®

The Sharpsafe® NPSD™ Range



The Near Patient Sharps Disposal System™
The safer system of sharps segregation

frontiermedical
group

Sharpsafe is a registered trade mark of Frontier Plastics Limited.
© June 2011 Frontier Medical Products Limited. All rights reserved. GP0059_A21_09

Frontier Medical Products
Newbridge Road Industrial Estate
Blackwood, South Wales NP12 2YN
Tel: +44(0)1495 235800
Fax: +44(0)1495 235808
Email: info@frontiermedical.eu

www.frontiermedical.eu

Reducing Sharps Waste;
The most cost effective and environmentally friendly way of
dealing with waste is to reduce the amount of waste produced

NetRegs2010

www.frontiermedical.eu

Sharpsafe®

The complete range of sharps containers



Sharps Containers		0.2L	0.3L	0.45L	0.6L	1L	2L	3L	4L	7L	9L	13L	24L	30L	Quiver 7.5L	25L XL	Clinisafe 8L	Clinisafe 11L
INFECTIOUS SHARPS FOR INCINERATION	Cat No. NHS No.	41731430	41751430 FSL359	41701430 FSL360	41502430 FSL336	41602430 FSL078	41405430 FSL426	41455430 FSL276	41005430 FSL030	41105430 FSL031	41172430 FSL290	41152430 FSL291	41202430 FSL308	41802430 FSL066	41902430 FSL080	41952430	41912430 FSL919	41922430
INFECTIOUS SHARPS FOR INCINERATION OR ALTERNATIVE TREATMENT	Cat No. NHS No.				41502410 FSL350	41602410 FSL337	41405410 FSL159	41455410 FSL351	41005410 FSL338	41105410 FSL353	41172410 FSL341	41152410 FSL100	41202410 FSL342	41802410 FSL354				
CYTO SHARPS FOR INCINERATION	Cat No. NHS No.				41502420 FSL335	41602420 FSL321	41405420 FSL322	41455420 FSL277	41005420 FSL326	41105420 FSL283	41172420 FSL327	41152420 FSL329	41202420 FSL334	41802420 FSL279				

FSL codes apply to NHS Supply Chain references only.

Leakproof Containers		6L	10L	19L
INFECTIOUS SHARPS FOR INCINERATION	Cat No.	41652430	41142430	41192430
CYTO SHARPS FOR INCINERATION	Cat No.	41652420	41142420	41192420
LIQUID PHARMA WASTE FOR INCINERATION	Cat No.	41652440	41142440	41192440


The Sharpsafe® Leakproof range is now available in 6 Litre, 10 Litre and 19 Litre containers which comply with both UN1851 and UN3291.

Solid Pharmaceutical Waste Containers		7L	30L
SOLID PHARMA WASTE FOR INCINERATION	Cat No. NHS No.	41105440 FSL284	41802440 FSL278

The Sharpsafe® Solid Pharmaceutical Waste container has been designed for POM's and other solid pharma waste in compliance with UN3249.

Specimen Transport Bags		STB1 6L	STB2 19L
COMMUNITY SPECIMEN TRANSPORT BAG	Cat No.	50010000	50020000

The Community Specimen Transport Bag (STB) has been designed for community healthcare workers who need to transport diagnostic specimens whilst adhering to UN3373.

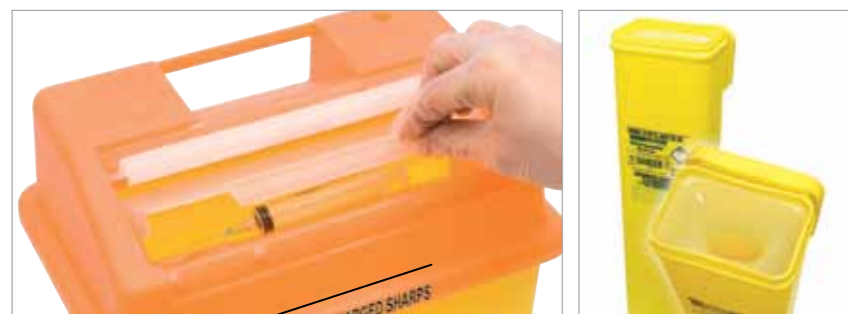

Protected Access Containers		9L	13L	24L	30L	Quiver 7.5L	25L XL
INFECTIOUS SHARPS FOR INCINERATION	Cat No. NHS No.	41171430 FSL292	41151430 FSL361	41201430 FSL363	41801430 FSL285	41901430	41951430
INFECTIOUS SHARPS FOR INCINERATION OR ALTERNATIVE TREATMENT	Cat No. NHS No.	41171410 FSL037	41151410 FSL362	41201410 FSL364	41801410 FSL355		
CYTO SHARPS FOR INCINERATION	Cat No. NHS No.	41171420 FSL328	41151420 FSL366	41201420 FSL365	41801420 FSL356		

Designed for use in public access areas such as A&E and out-patients. The Protected Access range allows easy disposal of sharps but provides added security for both user and the general public.

The Protected Access range also discourages the use of the sharps container for non sharps waste, avoiding unnecessary waste disposal costs.

NPSD™ - Helping to reduce needlestick injuries, reduce waste and reduce costs

Sharpsafe® SMS, Sharps Management System is a sharps containment and management system which incorporates Sharpsafe® Near Patient Sharps Disposal™ (NPSD™) and Aseptic Non Touch Technique (ANTT®).

NPSD™ Tray
Cat No. 43000000



What is NPSD™?

Near Patient Sharps Disposal™ (NPSD™) is a strategy for reducing needlestick injuries which features sharps containment products designed specifically to be taken to the patient. Sharps are disposed of safely at the 'point of use' thus eliminating many of the opportunities for needlestick injuries.

NPSD™ + ANTT : Aseptic Non Touch Technique

Aseptic Non Touch Technique (ANTT) is a DOH and EPIC2 endorsed 'best practice' aseptic technique developed to standardise practice, improve clinical standards and further protect patients and healthcare providers. The Sharpsafe® system incorporates trolleys and trays designed for use with ANTT.



Pedestal Trolley

Cat No. 44800000




Either or both




Sharpsafe® Trolley

Cat No. 44000000

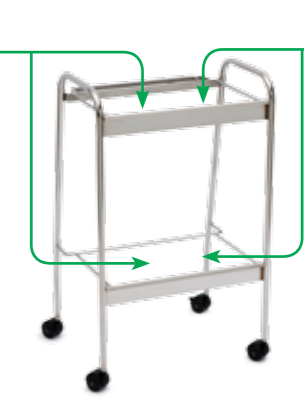

Trolley & Tray System

The Sharpsafe® Trolley can be configured with a variety of Sharpsafe containers and/or a variety of trays.

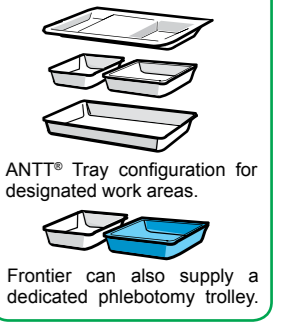

Clinical Waste Containers

30 Litre	Cat No. 41850000	60 Litre	Cat No. 41860000
----------	------------------	----------	------------------



Product Stands

Community Stands	Desktop Stands
------------------	----------------



Universal Bracket

Trolley Bracket
Cat No. 48000000



Medi-rail Bracket

Cat No. 48000001



Wall Bracket

Cat No. 48000002



Combination Pack Cat No. 48010000

Mounting Bracket
Cat No. 48000003

